

NuFlex Work-Surface-Mounted System

NuFlex adjustable shelving over Hanson fixed cabinets

- Combines the advantages of a fully integrated, pre-plumbed, pre-wired reagent rack system and a traditional floor-mounted laboratory casework format.
- Removable, adjustable reagent rack shelving provides more flexibility than standard fixed shelving.
- Utilities can be provided from overhead or under the countertop.
- The option to combine the advantages of a NuFlex free-standing floor-mounted system with standard fixed casework by utilizing the NuFlex work surface-mounted reagent rack system.

More NuFlex Benefits:

- **Reconfigure your lab as your needs change.** NuFlex is macro-adaptable and micro-adjustable, to easily accommodate new equipment, science, research, projects, processes, or researchers.
- **Does not lock in a fixed floor plan,** and returns control of reconfiguration to the system owner.
- **Optional bar-code inventory control system.** Makes it easy to move, change, and re-use NuFlex components while keeping track of all pieces.
- **Grow now, grow later.** We've been providing quality steel lab furniture to the science industry for over 30 years. We're with our customers for the long haul.
- Casework and shelving finishes available in **steel, wood veneer, plastic laminate, and phenolic.** Ask us about flexible finishes and custom colors.
- Modular NuFlex pieces can be combined with **fixed Hanson cabinets,** chem-resin work surfaces, fixtures and accessories.
- **Custom components** are also available — just ask.




NuFlex
THE SMART CHOICE FOR CHANGING LABS



NuFlex Flexible Furniture

- Flexible starter-adder system
- Exclusively Hanson: Unique sliding track for **Infinite, Anytime** adjustment of countertops and shelving — see more inside!
- Integrated data, plumbing, power, and fixtures dramatically reduce project costs
- Change it, add to it, move it, and take it with you
- Less costly per linear foot than other systems

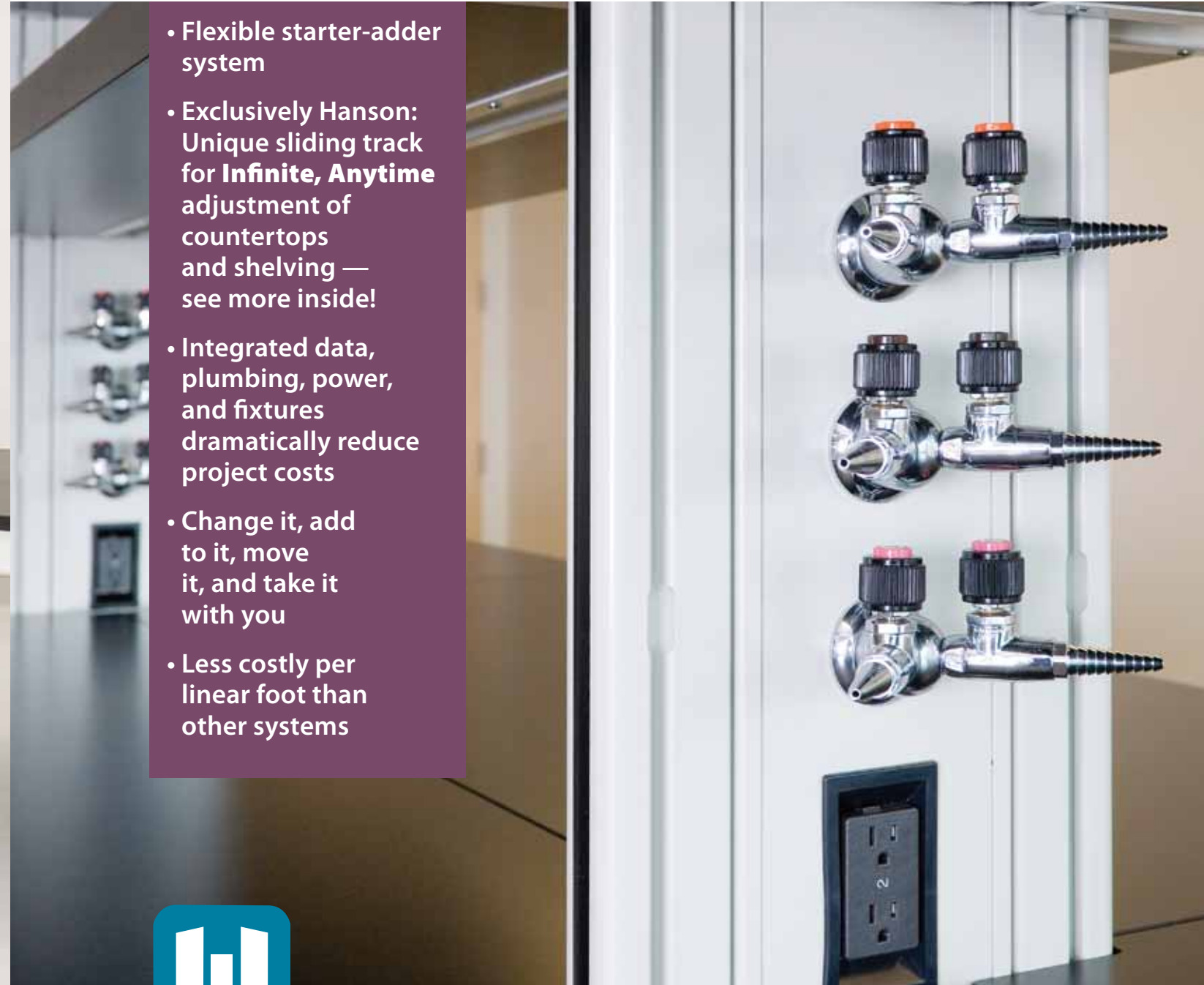

There's only one NuFlex

Because it's the only Flex you need

NuFlex Free-Standing Floor-Mounted System

Hanson NuFlex offers the ultimate combination of flexibility, strength, quality and features — for any lab, any application, now and in the future.

- Configure and reconfigure entire labs from the full line of NuFlex products.
- Reconfigure NuFlex in an existing lab, or take it with you to a new building.
- Taller, 84" uprights allow room for three shelves above the work surface, vs. standard two. Install countertops at any height from 26" to 46".
- Innovative design of upright supports offers huge advantages: it's cleaner and less susceptible to contamination. It's infinitely adjustable, so you can position your countertop and shelving anywhere. It looks neater and more attractive, too.
- Option for hydraulically adjustable countertop.

Exclusively Hanson

Hanson is the **only** lab furniture company to offer these unique features:


- **Unique sliding track** gives allows **Infinite, Anytime** adjustment of countertops and shelving. **Infinite:** loosen the retainer bolts, slide the surface to any position you like, and tighten the retainer bolts to lock in place. You're not stuck with the predetermined intervals of a ratchet system. **Anytime:** because you can change shelf or countertop position without dismantling the casework frame, it's fast and easy to make changes on the fly.
- Sliding track is **cleaner and less susceptible to contamination** than toothed ratchet system.
- **Add shelving below the work surface**, for maximum storage efficiency.
- You get a single point of connection for four 20-amp circuits. This **4-in-1 connection** is neater, cleaner, and more convenient.
- Position **electrical raceways** above or below the work surface.
- **One-inch pass-through** behind shelving and work surface allows neat, easy cabling.
- Three choices of shelving depth provide just the right usable space for your application. For the first time, you can select the shelving you need.
- Option to control work-surface height **manually or electrically**.



Call the Hanson Lab Planning Specialists

805-498-3121 or visit www.hansonlab.com/nuflex

for the full NuFlex product lines, sample configurations, and photos of NuFlex labs

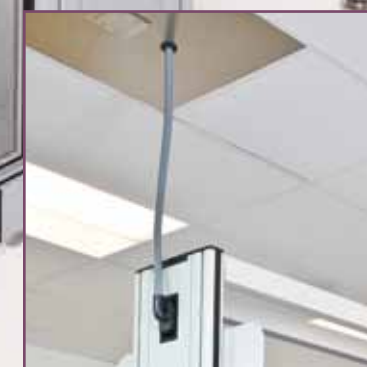

4-in-1 electrical connection



Specialty gas service connections



NuFlex offers multiple configurations

NuFlex
THE SMART CHOICE FOR CHANGING LABS



Hanson